
KiCA Evo-mini

Instructions

EN / CN / DE / FR / ES / IT /
KR / ID / PL

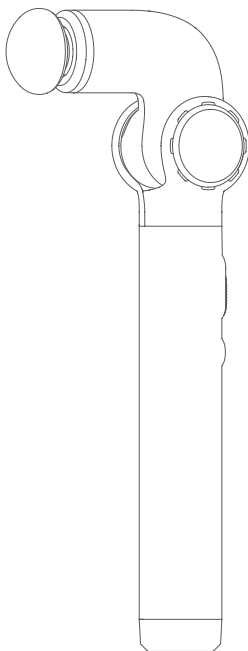- Warning: Please read the user manual carefully before use. Improper use of this product may cause damage and pose a risk to your safety.

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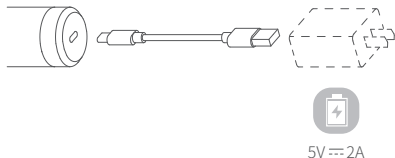
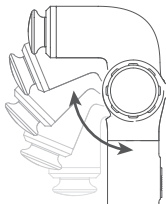
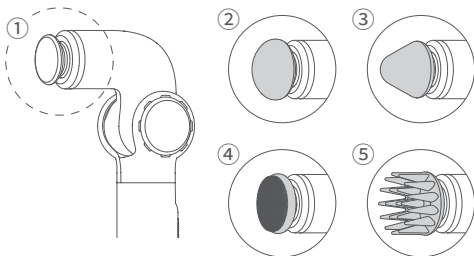
Flat head: Massage flat major muscle groups on the back.

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Spherical head: Relax soft tissues and sensitive muscles.

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Cone head: Percuss deep tissues like meridians and acupoint, etc.

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Gentle head: Enhance essential oil absorption.

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Hedgehog massage head: Mainly for the scalp, hairline, and neck care.

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Bouton d'alimentation

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
Tombol daya

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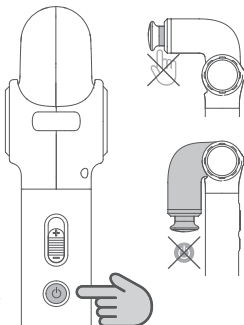Do not touch the middle and the back of the two massage heads while the product is running to avoid injuries. Do not power on the device when it is folded.

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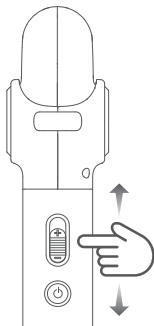
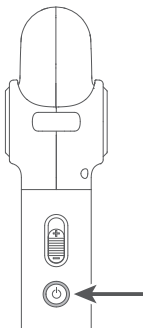
7

Intelligent protection

Mode	Description	Indicator
Continuous working protection	After working for 10 minutes continuously, the product will enter standby mode automatically (Single tap to wake the machine)	Green light flashes slowly
Locked motor protection	After the rotation is blocked, the machine will power off automatically	
High-temperature protection	The machine will power off automatically when the temperature is too high	Red and green lights flash alternately

Mode	Description	Indicator
Standby mode	If there is no operation for 3 minutes, the machine will power off automatically	Turn off after green light flashes slowly
Usage is restricted during charging	When connected to an external power source, the device stops running	Constant red light
Auto power off protection	—	Turn off after red light flashes quickly

Please do not press the massage head too hard to block the motor from rotating during use, as this can cause additional power consumption and cause the device to shut down due to low battery power, and can also reduce the service life of the product.

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【EN】

Contraindications:

1. Cancers
2. Metal bone anastomosis (do not use the device directly in the fracture area)
3. Osteoporosis
4. Vein thrombosis and pulsation
5. Hypertension
6. Pregnant woman
7. Varicose veins
8. Acute inflammation
9. Do not massage the bony area
10. Do not massage surgical site
11. Do not massage the cervical vertebrae and spine
12. Kidney stones (do not use the device near the kidneys)
13. Do not massage one episode for more than 5 minutes

Safety Instruction:

1. It is recommended to massage your body after wearing dry clothes or placing a clean and dry towel on massage part, then gently press and move the device.
2. High-speed and high-pressure massage can only be used to soft tissues without feeling pain or maladjustments. Do not massage head and any skeletal sites with high speed and high pressure.
3. Massage the same part of the body with high speed and high pressure for a long time will.
4. Please keep your fingers, hair or, any other body hair away from the extendable part of the device to prevent pinching and clipping.
5. Do not discard or misuse this device.
6. Do not modify this device by yourself without permission.
7. Do not operate or charge the device when nobody was around.
8. This device is for adults only.
9. People with contraindicated symptoms should use this device with doctor's permission.

Notice:

1. Make sure the massage head is not blocked by external force when the device is power on. If it is blocked, it will pause for a while and try to restart.
2. DO NOT contact the product with water or other liquid for the massage device is not waterproof.
3. DO NOT disassemble the product except marked detachable. It need send to after-sales or authorized service center to fix it if you accidentally disassemble it and cause abnormal work. The relevant costs are borne by users.
4. Prolonged continuous operation may cause the motor surface temperature to rise, please operate carefully.
5. DO NOT drop or strike the product. If the product is abnormal, contact After-sales support.

Storage and Maintenance:

1. Keep the product out of the reach of children and pets.
2. DO NOT leave the product near heat sources such as a furnace or heater. DO NOT leave the product inside of a vehicle on hot days.
3. Please store the product in a dry environment.
4. Do not overcharge or overdischarge the battery, otherwise it will damage the battery cells. It is recommended to charge at least once a month

when not in use.

5. Never use the product when the temperature is too high or too low.

Precautions when discarding the product:

The lithium ion battery used by the product is a recyclable and renewable valuable resource.

When discarding it, please have it handled by a professional agency that can recycle chargeable batteries. Before discarding the product, power it off and take the lithium ion battery out from it following the steps below. The battery is to be disposed of safely. Steps:

1. Disassemble the top housing with help of a professional tool.
2. Pull out the batter terminals from the PCBA.
3. Use a professional cutter to cut off the bottom housing and then take out the battery.

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